

Trace Elemental Species Separation and Detection



Contents

Introduction

The Ion Chromatography System	3
Speciation Analysis	4
Complete Inorganic Analysis Periodic Table of Total Inorganic Analysis	5
Arsenic IC-ICP-MS Speciation Analysis of As in Apple Juice Using the Thermo Scientific iCAP Q ICP-MS	6
ICP-ICP-MS Speciation Analysis of As in Organic Brown Rice Syrup (ORBS) Using the Thermo Scientific iCAP Q ICP-MS	7
Chromium Speciation Analysis of Cr(III) and Cr(VI) in Drinking Waters Using Anion-Exchange Chromatography Coupled to the Thermo Scientific iCAP Q ICP-MS	8
lodide Determination of Iodide and Iodate in Soy- and Milk-Based Infant Formulas	9
Mercury Combining the Synergies of Ion Chromatography and Inductively Coupled Plasma to Identify Mercury Contamination in Herbal Medicines	10
Total and Speciation Analysis of Mercury in Contact Lens Solutions by ICP-MS	11
Literature References	12

The Ion Chromatography System



Innovative IC Solutions for All Applications and Performance Needs

Thermo Scientific[™] Dionex[™] IC
Systems have led the analytical
instrument industry for over 30
years with solutions that
represent state-of-the art
technological advancements and
patented technologies. Our
products have evolved over many
generations, with each new
product providing enhanced
performance, greater reliability,
and easier operation.

Our High-Pressure Ion Chromatography HPIC™ Systems, including the Thermo Scientific Dionex ICS-5000⁺ system, are optimized for flexibility, modularity, and ease-of-use, combining the highest chromatographic resolution with convenience. Additionally, Capillary IC takes convenience to a whole new level.

The Thermo Scientific Dionex ICS-4000 is the world's first dedicated capillary high-pressure Reagent-Free™ (RFIC™) IC system, delivering high-pressure IC on demand with a system that is always ready for the next analysis.

Reagent-Free IC systems eliminate daily tasks of eluent and regenerant preparation in turn saving time, preventing errors, and increasing convenience. RFIC-EG systems use electrolytic technologies to generate eluent on demand from deionized water, and to suppress the eluent back to pure water to deliver unmatched sensitivity.

At the heart of our ion chromatography is a unique set of

column chemistries that provide high selectivities and efficiencies with excellent peak shape and resolution. Thermo Scientific™ Dionex[™] IonPac[™] chromatography columns address a variety of chromatographic separation modes including ion exchange, ion exclusion, reversed-phase ion pairing, and ion suppression. Our column chemistries are designed to solve specific applications, and we offer a variety of selectivities and capacities for simple and complex samples. Additionally, our Dionex IonPac column line is available in standard bore, microbore and capillary formats for the ultimate application flexibility. Learn more about our IC systems and consumables at www.thermoscientific.com/ dionex.



The complete Thermo Scientific Dionex IC Systems family

Speciation Analysis



Separate and Quantify Different Element Chemical Forms

The need to distinguish between chemical forms of an element has become critical for multiple industries, including the food, environmental, and pharmaceutical sectors. In the past, measuring the total amount of an element was sufficient. Unfortunately, the effects of an element extend far beyond its absolute amount. Different forms of an element can exhibit very different physiochemical properties, including varying toxicities. The process of separation and quantification of different chemical forms of an element, more specifically termed speciation analysis, can determine an element's various chemical forms. and thus deliver a better understanding of the environmental or health related impact associated with a particular sample. Speciation analysis can be split into two components: separation of individual ionic species by ion chromatography followed by trace elemental detection and quantification using ICP-MS. This combined method is termed Ion Chromatography-Inductively Coupled Plasma Mass Spectrometry (IC-ICP-MS).

ICP-MS is a multi-element spectrometry method to determine total elemental concentrations without bias towards metal species. The technique provides rapid and robust total element concentration determinations in various types of samples with high sensitivity (sub-part per thousand [ppt] detection). This technique leverages the combination of an ICP source with a mass spectrometer. The ICP is a high temperature source that decomposes and atomizes molecules, then ionizes the atoms. The mass spectrometer separates and detects these ions.

The process of ion chromatography allows ions and polar molecules to be separated on the basis of their charge, size, and polarizability. This specific method determines ionic species mainly with a conductivity detector, but can also be used with other types of detectors. With its metal-free fluidic flow path, ion chromatography is ideally suited to elemental speciation analysis. Furthermore, the system can analyze a range of compounds, from anionic and

cationic contaminants to disinfection by-products. These are all important indicators of quality in environmental waters, the pharmaceutical industry, and food applications; many of which are toxic and need to be regulated.

The Thermo Scientific $^{\scriptscriptstyle{\mathsf{T}}}$ iCAP $Q^{\scriptscriptstyle{\mathsf{T}}}$ ICP-MS system represents a unique platform to determine the total elemental concentration of a sample. The iCAP Q allows for high sensitivity that can provide single figure ppt detection limits for many elements. As a result, a full mass range analysis can be carried out for routine samples. Additionally, the iCap Q series houses a proprietary collision/ reaction cell with low mass cut-off so unwanted species do not pass through the quadrupole mass filter. When coupled with Dionex IC systems, these techniques successfully provide a complete picture for analyzing both total elemental concentration as well as chemical form of the element of interest.

Learn more about speciation analysis at www.thermoscientific.com/speciation.



The Thermo Scientific iCAP Q with the Thermo Scientific Dionex ICS-5000+ IC System

Complete Inorganic Analysis













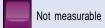

Periodic Table of Total Inorganic Analysis

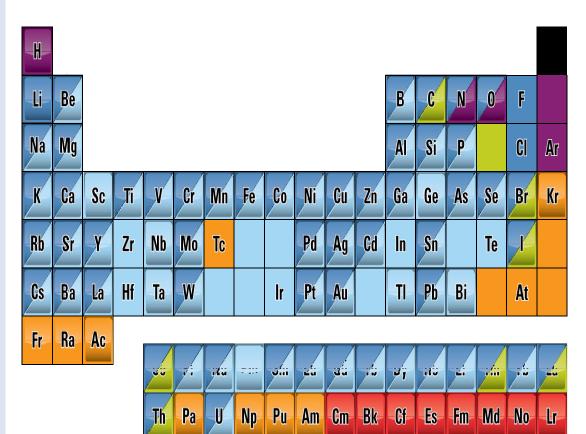

Análisis de especiación IC-ICP-MS de As en jugo de manzana utilizando el Thermo Scientific iCAP Q ICP-MS



Introducción

Informes recientes de los medios de comunicación en los Estados Unidos han afirmado que algunos jugos de manzana pueden contener altos niveles de arsénico. Sin embargo, en estos estudios de determinación, sólo se evaluó la concentración total de arsénico; no se llevó a cabo ninguna investigación detallada de la forma química del elemento

Esta es una distinción importante ya que las formas inorgánicas de arsénico (Como As(III) y Como As (V)) son altamente tóxicos, mientras que la forma orgánica (por ejemplo, arsenobetaína) no es considerado Por ser nocivo. Los niveles encontrados de arsénico total enel jugo de manzana son más bajos que el límite máximo de la US EPA (MCL) de 10 ng/g así que el jugo de manzana es considerado seguro y actualmente no está regulado.

Sin embargo, la FDA está revisando los datos actualmente y, puede reducir sus pautas actuales.

Download the Full Version of Application Note 43099

Arsenic Speciation in Apple Juice

Equipment

- Sistema de Dionex ICS-5000 IC *
- Columna Dionex IonPac AS7,2 mm i.d. × 250 mm
- iCAP Q ICP-MS
- *Dionex ICS-5000* Sistema HPIC puede ser usado para obtener resultados equivalentes

iCAP Q Operating Parameters		
Poder:	1550 W	
Gas Nebulizador	0.80 L/min	
Injector:	2 mm I.D., quartz	
Interface	Ni sampler and skimmer	
QCell He Flujo d	Gas: 4.8 mL/min	
QCell KED:	2 V	
75As Time: 10	0 ms	

Dionex ICS-5000 IC Operating Parameters		
Elución:	Gradiente	
Fase Movil:	A: 20 mmol/L Carbonato de Amonio B: 200 mmol/L Carbonato de Amonio	
Volume de Inyec	ción: 20 µL	
Duración:	15 min	

Analysis

IC-ICP-MS

Results

Ver tablas.

As species concentrations, method detection limits (MDLs) and total arsenic concentrations in two of the apple juice samples analyzed

	AsB	DMA	As³+	AsC	ММА	As5+	Sum of Total Species	Total As
MDL	0.002	0.004	0.005	0.004	0.011	0.001	-	0.005
Jugo 3	ND	ND	0.5 ± 0.01	ND	ND	0.7 ± 0.01	1.2	1.7 ± 0.05
Jugo 4	ND	0.4 ± 0.05	0.3 ± 0.01	ND	0.1 ± 0.05	0.7 ± 0.01	1.5	1.8 ± 0.05

All concentrations have units of ng/g. ND indicates not detected.

Spike recovery for six arsenic species in apple juice

Especies	Esperado (ng/g)	Encontrado (ng/g)	Recupración (%)
AsB	2.19	2.27	104
DMA	1.40	1.15	82
As (III)	1.35	1.38	102
AsC	1.94	1.87	94
MMA	1.09	1.13	104
As (V)	1.10	1.07	98

IC-ICP-MS Speciation Analysis of As in Organic Brown Rice Syrup (ORBS) using the Thermo Scientific iCAP Q ICP-MS



Introduction

Media reports and scientific publications on the determination of arsenic (As) in foodstuffs have sparked renewed interest from consumer groups and politicians leading to responses from national regulatory bodies. Based on recent reports, the FDA began carrying out a study on As in rice and rice products, including organic brown rice syrup, OBRS, an ingredient in a variety of organic foods. In this study, OBRS samples were analyzed for their total arsenic content by ICP-MS and then subsequently by IC-ICP-MS to determine the concentration of six arsenic species: the two toxic inorganic species As (III) and As (V), and four organic species that are considered harmless.

Download the Full Version of Application Note 43126

Arsenic Speciation in Organic Brown Rice Syrup

Equipment

- Dionex ICS-5000*
- Dionex IonPac AS7 column, 2 mm i.d. × 250 mm
- iCAP Q ICP-MS
- *Dionex ICS-5000* HPIC system can be used for equivalent results

iCAP Q Operating Parameters**		
1550 W		
0.80 L/min		
2 mm I.D., quartz		
Ni sampler and skimmer		
4.8 mL/min		
2 V		
100 ms		

Total As concentrations of >100 ng/g were found in three analyzed samples.

	⁷⁵ As (ng/g)
OBRS Sample #1	118 ± 7
OBRS Sample #2	136 ± 7
OBRS Sample #3	107 ± 11

Dionex ICS-5000 IC System Operating Parameters**		
Elution:	Gradient	
Mobile Phase:	A: 20 mmol/L Ammonium carbonate B: 200 mmol/L Ammonium carbonate	
Injection Volume:	20 μL	
Duration:	15 min	


**Operating parameters used in analysis based on configurations in AN 43126

Analysis

IC-ICP-MS

Results

See table and figure below.

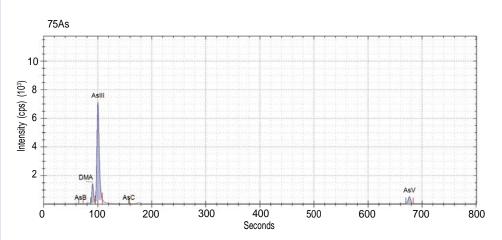

IC-ICP-MS chromatogram of arsenic species found in ORBS sample. As (III) was the most abundant species detected.

Speciation Analysis of Cr (III) and Cr (VI) in Drinking Waters Using Anion Exchange Chromatography Coupled to the Thermo Scientific iCAP Q ICP-MS



Introduction

Both the United States EPA and the European Union have specified maximum admissible chromium concentrations in their respective drinking water directives. As with many other trace elements, chromium (Cr) is typically found in more than one chemical form, each of which with different chemical properties and behavior, such as bioavailability and toxicity. For chromium, Cr (III) is essential to human beings and involved in different processes in the body while Cr (VI) is highly toxic. Total Cr content therefore in, for example, a drinking water sample does not provide sufficient information to evaluate potential hazards to populations exposed to it. In order to provide this critical information a supporting speciation analysis is required to determine the amounts of the different Cr species in the sample.

Download the Full Version of Application Note 43098

Chromium Speciation in Drinking Water

Equipment

- Dionex ICS-5000 IC system*
- Dionex IonPac AG7 column, 2 mm I.D. × 50 mm
- iCAP Qc ICP-MS
- *Dionex ICS-5000+ HPIC system can be used for equivalent results

iCAP Q Opera	ting Parameters
Forward Power:	1550 W
Nebulizer Gas:	0.80 L/min
Injector:	2 mm I.D.
QCell He Gas Flo	w: 4.8 mL/min
QCell KED:	2 V
Dwell Time:	100 ms

Dionex ICS-5000 IC System Operating Parameters


Elution:	Isocratic
Mobile Phase:	0.4 mol/L HNO ₃
Flow Rate:	400 μL/min
Injection Volume:	20 μL
Duration:	150 s

Analysis

IC-ICP-MS

Results

See figures and table below.

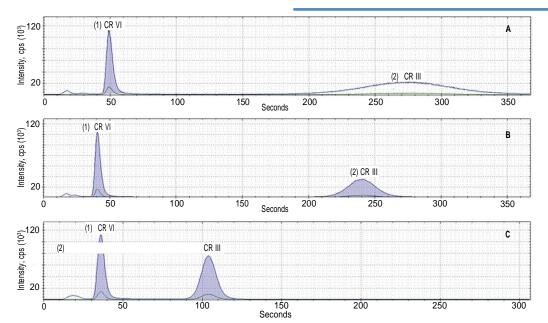

Cr(III) and Cr(VI) chromatograms obtained using 0.2 (A), 0.3 (B) and 0.4 (C) mol/L nitric acid as mobile phase. Note that the x-axis in (C) has been shortened to 300 s.

Recovery of Cr (VI) and (III) species from drinking water samples

Conc. spiked (ng/g)	Cr	(VI)	Cr	(III)
Conc. spiked (ng/g)	Found (ng/g)	Recovery (%)	Found (ng/g)	Recovery (%)
2.34 of each	2.31 ± 0.01	99 ± 1	2.35 ± 0.02	100 ± 1
6.03 Cr (VI); 1.90 Cr (III)	6.01 ± 0.02	100 ± 1	2.00 ± 0.01	105 ± 1
1.87 Cr (VI); 6.20 Cr (III)	1.85 ± 0.01	99 ± 1	6.15 ± 0.03	99 ± 1

Determination of Iodide and Iodate in Soy- and Milk-Based Infant Formulas



Introduction

Accurate measurement of iodine in food matrices requires a robust iodine extraction method and a sensitive analytical method for iodine quantification. This application note includes the acetic acid digestion method for iodide extraction, coupled with an IC-Pulsed Amperometric Detection (PAD) method for iodide detection first developed in an archived version of this application note. The IC method coupled with electrochemical detection allows for selective and sensitive determination of iodide in complex matrices. The acid digestion procedure to extract iodide was optimized for milk- and soy-based infant formulas. In addition, sample preparation conditions to convert iodate to iodide for determining total iodine (i.e., iodide and iodate) are presented.

Download the Full Version of Application Note 37

lodide and lodate Speciation in Infant Formula

Equipment

• Dionex ICS-5000 IC system* including: Gradient pump

DC Detector Chromatography Compartment ED Electrochemical Detector without cell (P/N 079830)

ED Electrode, Ag, with gasket and polishing kit (P/N 079856)

Ag/AgCl reference electrode (P/N 061879) Dionex AS or AS-AP Autosampler

- Dionex[™] IonPac AG11 guard, 4 × 50 mm (P/N 44078)
- Dionex IonPac AS11 analytical, 4 × 250 mm (P/N 44076)
- Thermo Scientific[™] Dionex[™] OnGuard[™] II RP Cartridges, 2.5 cc (P/N 057084)
- EO Eluent Organizer with two 2 L plastic bottles and pressure regulator
- Vial Kit, 0.3 mL polyprop with caps and septa (P/N 055428)
- Micro Tubes 1.5 mL, type D, without skirted base, screw cap assembled, sterile (Sarstedt[™] P/N 72.692.005 or equivalent)
- Thermo Scientific™ Nalgene™ narrow-mouth bottle, HDPE/PP, 1000 mL (P/N 2002-0032)
- Nalgene polystyrene lab filter unit, 500 mL upper capacity with 1000 mL receiver capacity 0.2 micron, 75 mm membrane diameter (P/N 154-0020)

Dionex ICS-5000 HPIC system can be used for equivalent results

Reagents and Standards

- Deionized (DI) water, type I reagent grade, $18\ M\Omega\text{-cm resistivity or better filtered through a} \\ 0.2\ \mu\text{m filter immediately before use}$
- · Nitric acid
- Sodium iodide
- Sodium iodate
- Ascorbic acid
- · Acetic acid


Dionex ICS-5000 IC System Operating Conditions		
Flow Rate:	1.5 mL/min	
Injection Volume:	100 μL	
Column Temp:	30°C	
Backpressure:	1000 psi	
Flush Volume:	1000 μL	
Detection:	PAD	
Cell Temp:	30°C	
Background:	2–10 nC	
Working Electrode:	Silver working electrode	
Reference Electro	de: Mode: Ag/AgCL Noise: 3–5 pC	


Analysis

IC-PAD

Results

See chromatograms below.

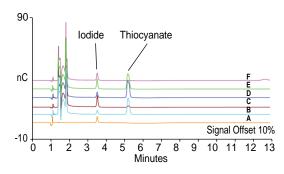

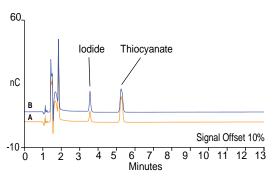

Left: Determination of iodide in DI water (A), milk-based infant formulas (B–E), and soy-based infant formula (F). Right: Chromatogram of iodide in infant formula (A) and infant formula spiked with iodate (B).

Combining the Synergies of Ion Chromatography and Inductively Coupled Plasma to Identify Mercury Contamination in Herbal Medicine



Introduction

Both mercury and lead are neurologic toxins and bioaccumulators, targeting brain, and other organs which can cause birth defects and sometimes death. Therefore, not only must the patients be rapidly diagnosed and treated to minimize the damage to their health but the source of their illness must be also rapidly identified to prevent other future cases. This application note demonstrates the advantages of using ion chromatography with inductively coupled plasma mass spectrometry for mercury speciation.

Download the Full Version of Application Note 43130

Mercury Speciation in Herbal Medicines

Equipment

- Dionex ICS-1600 IC System including:
 Dionex AS-AP Autosampler (P/N 074921)
 VWD Variable Wavelength Detector 3400 (P/N 070221)
- Dionex IonPac CS5A mixed cation/anion exchange column set
- Inductively Coupled Plasma Mass Specrometer (ICP-MS)

Standards and Reagents

- Acetic acid
- · Sodium perchlorate
- Cystine

Analysis

IC-ICP-MS

Results

In this mercury poisoning cluster example, the IC analysis provided mercury speciation which defines to potential toxicity based on the toxicity of each species, whereas the ICP-MS analysis provided a fast multi-element screening. This revealed that mercury caused the clinical symptoms, and determined the total mercury contamination. Inorganic mercury determinations are shown in the table below.

Results of mercury determinations in contaminated herbal medicines

	Inorganic Mercury					
Sample	Measured (mg/L)	Calculated (wt %)	Calculated (µg/tablet)	Methylmercury / Ethylmercury		
Control	_	-	_	-		
1	3.69	0.0586	52	ND		
2	15.6	0.2265	387	ND		
3	4.82	0.0437	117	ND		
4	5.31	0.0473	123	ND		
5	11.3	0.0708	243	ND		

Total and Speciation Analysis of Mercury in Contact Lens Solutions by ICP-MS



Introduction

While there is continual awareness regarding exposure to mercury (Hg) sources in general and MeHg* in particular due to its presence in food samples such as fish, less interest is paid to the potential risk from ethylmercury (EtHg+ or EtHgX). One of the main reasons for this is the faster degradation and consequently excretion of EtHg+ in the human body that results in considerably lower chronic toxicity. There remains however potential sources where acute intake of EtHg+ can occur, for example as a consequence of exposure to thiomersal. Thiomersal is used as a

bactericide in multi-dose and in other health related products such as eye drops or contact lens solutions. The compound hydrolyzes in aqueous solution to form EtHg⁺ and thiosalycilate which is an effective bactericide.

Download the Full Version of Application Note 43141

Mercury Speciation in Contact Lens Solutions

Equipment

- Dionex ICS-5000 IC System*
- Dionex IonPac CS5A column, 2 mm I.D. × 250 mm
- iCAP Qc ICP-MS
- *Dionex ICS-5000* HPIC System can be used for equivalent results

iCAP Q Operating Parameters				
Forward Power:	1550 W			
Nebulizer Gas:	1.05 L/min			
Injector:	2 mm I.D., quartz			
Interface:	Ni sampler and skimmer			
Dwell Time:	10 ms, 100 ms for speciation analysis			
Analysis Mode:	Standard (no cell gas)			

Dionex ICS-5000 IC System Operating Parameters				
Elution:	Isocratic			
Mobile Phase:	10 mmol/L NaClO ₄ 10 mmol/L acetic acid 10 mmol/L cystine			
Flow Rate:	0.5 mL/min			
Injection Volume:	20 μL			
Duration:	5 min			

Analysis


IC-ICP-MS

Results

See table and chromatogram below.

Spike recovery of thiomersal in contact lens solution

Sample #	Amount Spiked (mg/kg)	Amount Recovered (mg/kg)	Spike Recovery (%)
1	10.2	10.9 ± 0.04	108
2	18.1	18.8 ± 0.07	104

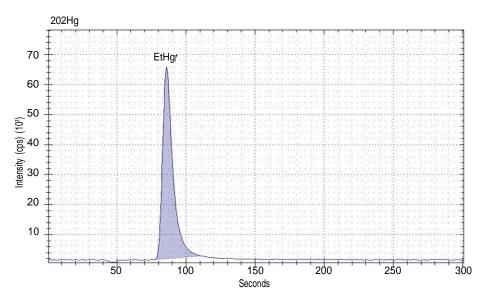

Chromatographic separation of EtHg+ derived from thimerosal hydrolysis

Literature References

References for Application Note 43099

- Letters from the FDA to the Dr. Oz Show Regarding Apple Juice and Arsenic (http://www.fda.gov/Food/ ResourcesForYou/Consumers/ucm271746.htm)
- FDA arsenic in apple juice resources site: http://www.fda.gov/Food/ResourcesForYou/Consumers/ucm271595.htm
- FDA arsenic in apple juice results: http://www.fda.gov/ Food/FoodSafety/FoodContaminantsAdulteration/Metals/ ucm272705.htm
- Dionex homepage (http://www.dionex.com/en-us/products/ columns/ic-rfic/specialty-packed/ionpac-as7/lp-73274.html)

References for Application Note 43126

- Arsenic, Organic Foods, and Brown Rice Syrup. Brian P. Jackson¹, Vivien F. Taylor¹, Margaret R. Karagas², Tracy Punshon³, Kathryn L. Cottingham³, ¹Trace Element Analysis Laboratory, Department of Earth Sciences, Dartmouth College, Hanover, NH, 03755, USA, ²Dartmouth Medical School, Section of Biostatistics and Epidemiology, Department of Community and Family Medicine, One Medical Center Drive, Lebanon, NH 03756, ³Department of Biological Sciences, Dartmouth College, Hanover, NH, 03755, USA. Environmental Health Perspectives, 2012, http://dx.doi.org/10.1289/ehp.1104619
- FDA Statement on Arsenic in Brown Rice Syrup, February 2012: http://www.fda.gov/Food/FoodSafety/FoodContaminantsAdulteration/Metals/ucm292531.htm
- IC-ICP-MS speciation analysis of As in apple juice using the Thermo Scientific iCAP Q ICP-MS, Application Note AN43099.
- Quantitative chemical extraction for arsenic speciation in rice grains. Huang et al., *J. Anal. At. Spectrom.* 2010, 25, 800–802. DOI: 10.1039/C002306.

References for Application Note 43098

- Séby, F., Charles, S., Gagean, M., Garraud, H., Donard, O. F. X., *J. Anal. At. Spectrom.*, **2003**, *18*, 1386–1390
- Xing, L., Beauchemin, D., J. Anal. At. Spectrom. 2010, 25, 1046–1055
- Dionex homepage (http://www.dionex.com/en-us/products/ columns/ic-rfic/specialty-packed/ionpac-as7/lp-73274.html)

References for Application Note 37

- Haldiman, M.; Alt, A.; Blanc, K.; Blondeau, K. Iodine Content of Food Groups. *J. Food Comp. Anal.* **2005**, *18*, 461–471.
- Dietary Reference Intakes for Vitamin A, Vitamin K, Arsenic, Boron, Chromium, Copper, Iodine, Iron, Manganese, Molybdenum, Nickel, Silicon, Vanadium, and Zinc; A Report of the Food and Nutrition Board, Institute of Medicine; National Academic Press, Washington, DC, 2001.
- Azizi, F.; Hedayati, M.; Rahmani, M.; Sheikloleslam, R.; Allahverdian, S.; Salarkia, N. Reappraisal of the Risk of Iodine-Induced Hyperthyroidism: An Epidemiological Population Survey.
 J. Endocrinol. Invest. 2005, 28, 23–29.
- Anke, M.; Groppel, B.; Müller, M.; Scholz, E.; Kramer, K. The Iodine Supply of Humans Depending on Site, Food Offer and Water Supply. Fresenius' J. Anal. Chem. 1995, 352, 97–101.
- Human Vitamin and Mineral Requirements; Report of a Joint FAO/WHO Expert Consultation: Bangkok, Thailand, 2002, Chapter 12. Iodine. www.fao.org/docrep/004/Y2809E/y2809e0i. htm#bm18 (accessed on Dec 7, 2012).
- Fisher, D.A. Upper Limit of Iodine in Infant Formulas. *J. Nutr.* **1989**, *119*, 1865–1868. http://jn.nutrition.org/content/119/12_Suppl/1865.full.pdf (accessed on Dec 7, 2012).
- Benofti, J.; Benotti, N. Protein-Bound Iodine, Total Iodine, and Butanol-Extractable Iodine by Partial Automation. *Clin. Chem.* 1963, 9 (4), 408–416.
- Mitsuhashi. T; Kaneda, Y. Gas Chromatographic Determination of Total Iodine in Foods. *J. AOAC Int.* **1990**, *73* (5), 790–792.
- Hammer, D.; Andrey, D. Comparison of Ion-Selective Electrode and Inductively Coupled Plasma-Mass Spectrometry to Determine Iodine in Milk-Based Nutritional Products. *J. AOAC Int.* 2008, 91 (6), 1397–1401.
- Crecelius, E.A. Determination of Total Iodine in Milk by X-ray Fluorescence Spectrometry and Iodide Electrode. *Anal. Chem.* 1975, 47, 2034–2035.
- Naozuka, J.; Messquita Silva da Veiga, A.M.; Oliveira, P.V.; de Oliveira, E. Determination of Chlorine, Bromine and Iodine in Milk Samples by ICP-OES. J. Anal. At. Spectrom. 2003, 18, 917–921.

Literature References

- Pacquette, L.H.; Levenson, A.M.; Thompson, J.J. Determination of Total Iodine in Infant Formula and Nutritional Products by Inductively Coupled Plasma/ Mass Spectrometry: Single-Laboratory Validation, J. AOAC Int. 2012, 95 (1), 169–176.
- Osterc, A.; Stos, K.; Stibilj. V. Investigation of Iodine in Infant Starting, Special and Follow-On Formulae. Food Control 2006, 17, 522–526.
- Niemann, R. A.; Anderson, D. L. Determination of Iodide and Thiocyanate in Powdered Milk and Infant Formula by On-Line Enrichment Ion Chromatography with Photodiode Array Detection. J. Chromatogr., A 2008, 1200, 193–197.
- Chadha, R.K.; Lawrence, J.F. Determination of Iodide in Dairy Products and Table Salt by Ion Chromatography with Electrochemical Detection. *J. Chromatogr.* 1990, 518, 268–272.
- Liang, L.; Cai, Y.; Mou, S.; Cheng, J. Comparisons of Disposable and Conventional Silver Working Electrode for the Determination of Iodide Using High-Performance Anion-Exchange Chromatography with Pulsed Amperometric Detection. *J. Chromatogr., A* 2005, 1085, 37–41.
- Dionex (now part of Thermo Scientific) Application Update 122:
 The Determination of Iodide in Brine, Sunnyvale, CA [Online]
 www.dionex.com/en-us/webdocs/4185-AU122_LPN034075-01.pdf
 (accessed Dec 10, 2012).
- Yang, S.X.; Fu, S.J.; Wang, M.L. Determination of Trace Iodine in Food and Biological Samples by Cathodic Stripping Voltammetry. *Anal. Chem.* 1991, 63, 2970–2973.
- Yebra, M.C.; Bollaín, M.H. A Simple Indirect Automatic Method to Determine Total Iodine in Milk Products by Flame Atomic Absorption Spectrometry. *Talanta* **2010**, *82* (2), 828–833.
- Official Methods of Analysis of AOAC International, 18th Ed., Official Method 992.24, *Iodine in Ready-to-Feed Milk-Based Infant Formula—Ion-Selective Electrode Method*, AOAC International, Gaithersburg, MD, 2007.
- Federal Ministry for Consumer Protection (1998) Amtliche Sammlung Untersuchungsverfahren nach § 35 LMBG 49.00-6, 1-4.
- European Committee for Standardization (2007) EVS-EN 15111:2007, Brussels, Belgium 1–12 [Online] www.evs.ee/preview/ evs-en-15111-2007-en. pdf (accessed Dec 10, 2012).
- International Organization for Standardization (2009) ISO 14378:2009, Geneva, Switzerland, 1–9.

- Knapp, G.; Maichin, B.; Fecher, P.; Hasse, S. Iodine Determination in Biological Materials Options for Sample Preparation and Final Determination. *Fresenius' J. Anal. Chem.* **1998**, *362*, 508–513.
- DRAFT AOAC Standard Method Performance Requirements (SMPR) 2012.XXX; Version 2; June 8, 2012: Determination of Iodine in Infant and Adult/Pediatric Nutritional Formula [Online] www.aoac.org/SPIFAN/Iodine_SMPR_v2.pdf (accessed Dec 10, 2012).
- Wang, L.M.; Lin, Y.K.; Wang, B.T.; Yan, Z.; Li, Y. Determination of Iodine Content in Infant Formula Cereals by Ion Chromatography. J. Instrumental Analysis 2011, 30 (1), 99–102.
- Sanchez, L.F.; Szpunar, J.; Speciation Analysis for Iodine in Milk by Size-Exclusion Chromatography with Inductively Coupled Plasma Mass Spectrometric Detection (SEC-ICP MS). *J. Anal. At. Spectrom.* **1999**, *14*, 1697–1702.

References for Application Note 43130

- Sarzanini, C.; Sacchero, G.; Aceto, M.; Aboilino, O.; Mentasti, Ed. Ion Chromatographic Separation and On-Line Cold Vapour Atomic Absorption Spectrometric Determination of Methylmercury, Ethylmercury and Inorganic Mercury, *Anal. Chim. Acta*, **1993**, *00*, 1–7.
- Dickson, H.R; Price, R. Application Note AN40992 Accurate Analysis of Low Levels of Mercury in Fish by Vapor Generation AA, Thermo Fisher Scientific, Cambridge, UK, 2010.
- Christison, T.; Hoefler, F.; Lopez, L. Application Note AN43130
 Combining the Synergies of Ion Chromatography and Inductively
 Coupled Plasma to Identify Mercury Contamination in Herbal
 Medicines, Thermo Fisher Scientific, Sunnyvale, CA, 2012.
- Dionex Application Note AN 131, Determination of Transition Metals at PPT Levels in High Purity Water and SC-2 (D-Clean) Baths, LPN 1058, Thermo Fisher Scientific, Sunnyvale, CA, 1998.

References for Application Note 43141

- Thimerosal in Vaccines Questions and Answers, US Food and Drug Administration; http://www.fda.gov/BiologicsBloodVaccines/ Vaccines/QuestionsaboutVaccines/UCM070430
- Statement on thiomersal, World Health Organisation; http://www.who.int/vaccine_safety/topics/thiomersal/ statement_jul2006/en/index.html
- Thimerosal in Vaccines, US Food and Drug Administration; http://www.fda.gov/BiologicsBloodVaccines/SafetyAvailability/ VaccineSafety/UCM096228

www.thermoscientific.com

©2013 Thermo Fisher Scientific Inc. All rights reserved. ISO is a trademark of the International Standards Organization. Sarstedt is a trademark of Sarstedt AG & Co.All other trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. This information is presented as an example of the capabilities of Thermo Fisher Scientific Inc. products. It is not intended to encourage use of these products in any manners that might infringe the intellectual property rights of others. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.



Thermo Fisher Scientific, Sunnyvale, CA USA is ISO 9001:2008 Certified.

Australia +61 3 9757 4486 Austria +43 1 333 50 34 0 Belgium +32 53 73 42 41 Brazil +55 11 3731 5140 China +852 2428 3282 Denmark +45 70 23 62 60 France +33 1 60 92 48 00 Germany +49 6126 991 0 India +91 22 6742 9494 Italy +39 02 51 62 1267 Japan +81 6 6885 1213 Korea +82 2 3420 8600 Netherlands +31 76 579 55 55 Singapore +65 6289 1190 Sweden +46 8 473 3380 Switzerland +41 62 205 9966 Taiwan +886 2 8751 6655 UK/Ireland +44 1442 233555 USA and Canada +847 295 7500



Part of Thermo Fisher Scientific